

## RAYCAST® HIGH PRECISION BASE PLATE-TO-COUCH TOP FIXATION DEVICES ADAPTER PLATES



Article No. : **35747**



Article No. : **35747/4**



Article No. : **35747/6**



Article No. : **35747/8**



Article No. : **32063**

### A. GENERAL PRODUCT INFORMATION

The products referred to in this document are devices that connect Orfit radiation therapy base plates to couch top fixation devices like the table fix and two or three-pin indexing bars.

This product may only be used in combination with immobilization masks produced by Orfit. Orfit prohibits the use of unauthorized third-party products in conjunction with its own products.


### B. PRODUCT RANGE

Art. No. 35747	Adapter plate for Table Fix – 12mm Base Plate (with screws) <i>Except 32113 &amp; 32117 Cyberknife base plates</i>
Art. No. 35747/4	Caudal Adapter plate for 2-pin Indexing Bar – 30mm Base Plate <i>For use only with two-pin bars</i>
Art. No. 35747/6	Adapter plate for 2-pin Indexing Bar – 12mm Base Plate <i>For use only with two-pin bars.</i>
Art. No. 35747/8	Caudal adapter plate for 2-pin Indexing Bar – 30 Base Plate (MRI Compatible) <i>For use with two-pin and three-pin indexing bar.</i>
Art. No. 32063	Adapter plate for 2-pin Indexing Bar – 12mm Base Plate (MRI Compatible) <i>For use with two-pin and three-pin indexing bars.</i>

### C. ASSEMBLY

#### Art. No. 35747

Place the Table Fix onto the treatment couch and then place the base plate – adapter plate assembly on top. Use the four M5 x 10 bolts and Allen wrench to fix the adapter plate to the Table Fix.

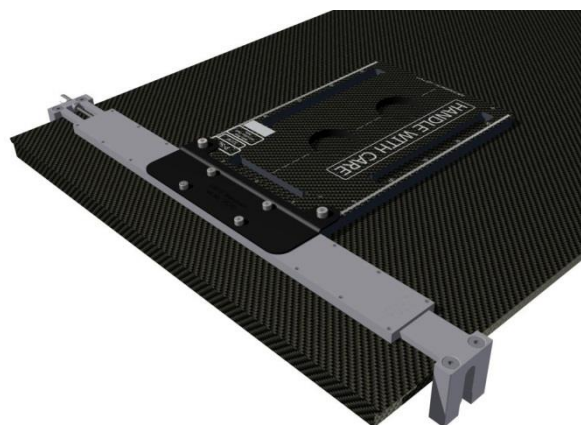

#### Parts list:

- 35747 – Adapter plate (Aluminium)
- 35744/7 – M5 x 10 bolt (4x)
- 35744/3 – Allen wrench 4 mm (1x)
- 35744/4 – M6 x 10 bolt (2x)
- 35744/5 – Allen wrench 5 mm (1x)

**Art. No. 35747/4**



*Note: The drawing shows a black adapter plate. In reality the adapter plate has a metallic grey colour.*

**Parts list:**

- 35747/4 – Adapter Plate (RVS)
- 35744/1 – M5 x 20 Screw (2x)
- 35744/3 – Allen Wrench 4 mm

Attach the adapter plate to the caudal end of the large base plates (AIO and IMRT long board) by means of the two M5 screws. This assembly can then be put over the pins of a two-pin lock bar.

Use together with other fixation devices like the AIO Loxon or the AIO fix as lateral fixators on the large base plates.

**Remark:**

Do not use this solution when the base plate is overhanging the couch top. It is recommended to use a table fix solution with L-shaped blocks to prevent tilting.

**Art. No. 35747/6**



*Note: The drawing shows a black adapter plate. In reality the adapter plate has a metallic grey colour.*

**Parts list:**

- 35747/6 – Adapter Plate (RVS)
- 35744/4 – M6 x 10 Screw (2x)
- 35744/5 – Allen Wrench 5 mm

Attach the adapter plate to the superior end of the base plate by means of the two M6 screws. This assembly can then be put over the pins of a two-pin lock bar.

Verify that the adapter plate is attached in the correct manner. The shape of the adapter plate should follow the shape of the base plate, where the adapter plate is bent.

**Art. No. 35747/8**



**Parts list:**

- 35747/8 – Adapter Plate (POM)
- 32066/2 – M5 x 20 Screw (2x)
- 35744/3 – Allen Wrench 4 mm


Attach the adapter plate to the caudal end of the large base plates (AIO and IMRT long board) by means of the two M5 screws. This assembly can then be put over the pins of a two or three-pin lock bar.

Use together with other fixation devices like the AIO Loxon or the AIO fix as lateral fixators on the large base plates.

**Remark:**

Do not use this solution when the base plate is overhanging the couch top. It is recommended to use a table fix solution with L-shaped blocks to prevent tilting.

**Art. No. 32063**

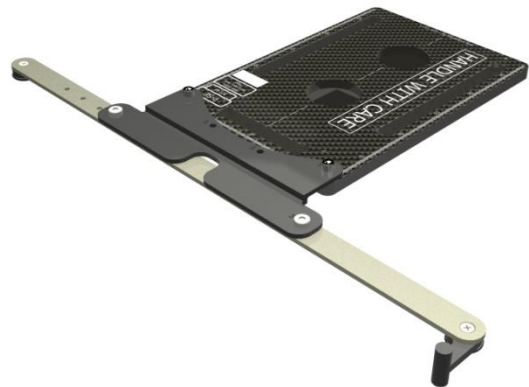

**Parts list:**

- 32063 – Adapter Plate (POM)
- 32063/1 – M6 x 12 Screw (2x)
- 35744/5 – Allen Wrench 5 mm

Attach the adapter plate to the superior end of the base plate by means of the two M6 screws. This assembly can then be put over the pins of a two or three-pin lock bar.

Verify that the adapter plate is attached in the correct manner. The shape of the adapter plate should follow the shape of the base plate, where the adapter plate is bent.

#### D. STORAGE

The base plate can be stored with adapter plate attached. If so, be careful during handling as the device is an extension to the base plate and may injure persons. Always store the item in such a way that the adapter plate does not take any weight of the base plate. Do not store anything on top of the adapter plate. A damaged adapter plate must be repaired or replaced!

Always store the product in a safe place to prevent it from getting damaged.  
Handle with care.

#### E. MAINTENANCE AND WASTE MANAGEMENT

This product can be cleaned and disinfected by means of soapy water or an isopropanol based disinfectant, applied with a soft cloth. If unsure about the cleaning fluid, do not use. **Never use aerosol sprays, corrosive cleaning agents, solvents or abrasive detergents.**

Periodic checks of the product should be done to insure the parts are not worn and require repair or replacement. **Do not attempt to make repairs yourself.** Contact your distributor if there are any questions or concerns.

The product can be disposed of with household waste.

#### F. ADDITIONAL INFORMATION

For additional information such as distributor contact information, product brochures, Safety Data Sheets and regulatory information, please visit our website [www.orfit.com](http://www.orfit.com).

#### Note:

The instructions were written in accordance with the European Directive 93/42/EEC for Medical Devices. It is prohibited to make alterations to this text without prior approval from Orfit Industries.

RAYCAST™ is a registered trademark of ORFIT INDUSTRIES N.V..



**ORFIT INDUSTRIES**  
Vosveld 9A | B-2110 Wijnegem | Belgium  
T (+32) (0)3 326 20 26  
[welcome@orfit.com](mailto:welcome@orfit.com)

[www.orfit.com](http://www.orfit.com)



Ref. No. 31937  
VERSION 12  
LAST UPDATE: 13/12/2017  
REVISION DATE: 13/12/2019