


THE AIO SOLUTION® 3.0 CARBON FIBRE BASE PLATE SBRT COMPATIBLE

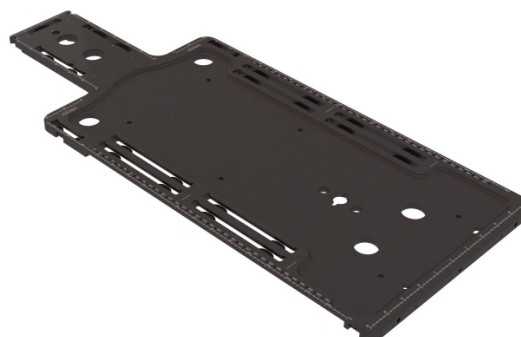

Article No. : **38013 & 38016**

A. GENERAL PRODUCT INFORMATION

This product is a medical device used for patient positioning and immobilisation in radiation therapy. The product can be used during both the simulation and treatment stage.

These products may only be used in combination with immobilisation masks produced by Orfit. Orfit prohibits the use of unauthorised third-party products in conjunction with its own products.

B. PRODUCT DESCRIPTION

The AIO 3.0 SBRT compatible base plate is the basic element of the AIO Solution 3.0. The base plate is used in combination with the AIO Solution 3.0 cushions and accessories, the AIO 3.0 SBRT bridges, the Raycast® High Precision Head Supports and Blocks & Wedges, and with the Efficast® Pre-cuts to form a reproducible patient positioning and immobilisation device in the field of radiotherapy. Information on these other parts and instructions on how to make the masks can be found in the respective 'instructions for use' and on www.orfit.com.

Pelvic mask slots on base plate 38013 and 38016 have a different position compared to other AIO 3.0 base plates. This makes base plate 38013 and 38016 not interchangeable with other base plates for use of lung, breast and pelvic masks.

SBRT pressure bridges (38218 & 38222) can be mounted on both edges of the AIO 3.0 SBRT compatible base plate. More information about the use and indexing of these pressure bridges can be found in the respective 'instructions for use'.

Hand holds (25000/7/A, 25000/7/B, 25000/69/A & 25000/69/B) can be mounted on both edges of the AIO 3.0 SBRT compatible base plate. More information about the use and indexing of these hand holds can be found in the respective 'instructions for use'.

Base plate 38016 has no leg separator insert which makes this base plate more homogeneous and low density in the pelvic area.

C. PRECAUTIONS FOR USE

This is a fragile product, please handle with care. Always place the base plate on a flat surface. Clear the treatment table of any debris before positioning the base plate. Always fix the base plate securely to the treatment couch before positioning a patient. More information on how to fix these base plates to a couch top can be found in the 'instructions for use' of the fixation devices.

The carbon fibre base plate is constructed to be light in weight and to have superb dosimetric properties. When handled roughly, it may get damaged and fibres may come off. When this happens, stop using the base plate to prevent fibres from getting in contact with the patient's or user's skin. Contact your distributor.

Do not lean or sit on the outer edge of the base plate. Do not apply point loads on the base plate as this might damage the base plate.

D. STORAGE

Always store the product in a safe place to prevent it from getting damaged or falling onto other objects. Take care not to damage the edges of the base plate when storing it in an upright position. Prevent hard objects from falling onto the base plate.

Store the base plate between +10°C (50°F) and 40°C (104°F).

E. PROPERTIES

38013 – AIO 3.0 – Base Plate CFL – SBRT compatible

Physical properties:

This plate is made of carbon fibre laminate.

Dimensions: L 1298 mm x W 530 mm x H 21 mm
L 51.10" x W 20.87" x H 0.83"

Weight: 3.88 kg



This product is MR unsafe and cannot be used in a MRI scanner.

Mechanical properties:


Maximum allowed patient weight: 200 kg

Maximum load on the head part when used in the maximum overhang position: 80 kg

Dosimetric properties:

Position	Attenuation (± 0.15 %) (6MV)	Attenuation (± 0.15 %) (15MV)	Skin build-up (± 0.1 mm) (15 MV)
	%	%	Mm H ₂ O equiv.
1	1.8	1.1	4.1
2	4.5	2.8	10.7
3	1.8	1.2	4.3
4	3.8	2.5	9.3
5	1.9	1.2	4.3
6	1.1	0.7	2.3
7	1.1	0.7	2.6
8	1.1	0.7	2.4
9	3.7	2.4	8.9
10	2.2	1.4	5.0

The drawing below shows the corresponding positions. The AIO 3.0 SBRT compatible base plate has areas where higher density materials have been used to reinforce parts like screw inserts and the leg separator insert. The base plate has low attenuation and skin build-up values in the treatment areas.

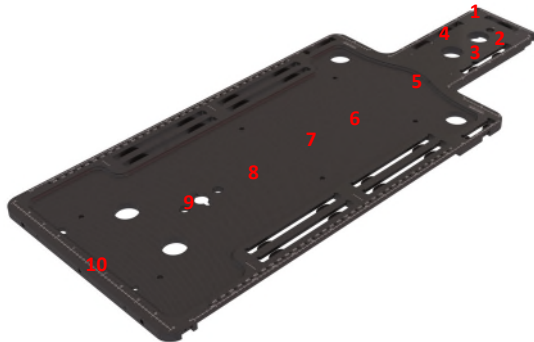

Note: Use these numbers as a guidance only. Perform the measurements again in your department according to applied treatment techniques and treatments systems.

38016 – AIO 3.0 – Base Plate CFL – SBRT compatible – Without Leg Separator

Physical properties:

This plate is made of carbon fibre laminate.

Dimensions: L 1298 mm x W 530 mm x H 21 mm
L 51.10" x W 20.87" x H 0.83"

Weight: 3.98 kg



This product is MR unsafe and cannot be used in a MRI scanner.


Mechanical properties:

Maximum allowed patient weight: 200 kg
Maximum load on the head part when used in the maximum overhang position: 80 kg

Dosimetric properties:

Values can differ from the 38013 – AIO 3.0 – Base Plate CFL – SBRT compatible. Perform the measurements again in your department according to applied treatment techniques and treatments systems.

The drawing below shows the corresponding positions. The AIO 3.0 SBRT compatible base plate has areas where higher density materials have been used to reinforce parts like screw inserts and the leg separator insert. The base plate has low attenuation and skin build-up values in the treatment areas.

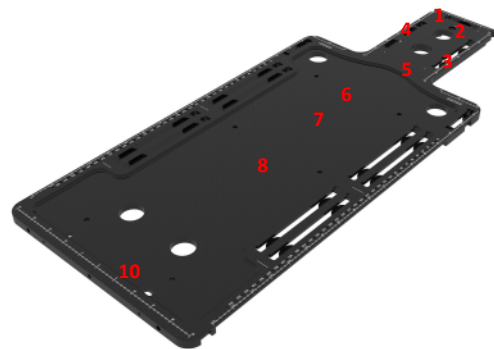

F. MAINTENANCE AND WASTE MANAGEMENT

This product can be cleaned and disinfected by means of an isopropanol based disinfectant, applied with a soft cloth. If unsure about the cleaning fluid, do not use. **Never use aerosol sprays, corrosive cleaning agents, solvents or abrasive detergents.** Further cleaning instructions can be found in the [Orfit Cleaning Guidelines](#).

Periodic checks of these products should be done to insure the parts are not worn and require repair or replacement. **Do not attempt to make repairs yourself.** Contact your distributor if there are any questions or concerns.

These products can be disposed of with household waste.

G. ADDITIONAL INFORMATION

For additional information such as distributor contact information, product brochures, Safety Data Sheets and regulatory information, please visit our website www.orfit.com.

Note:

THE AIO SOLUTION® is a registered trademark of Orfit Industries.



ORFIT INDUSTRIES
Vosveld 9A | B-2110 Wijnegem | Belgium
T (+32) (0)3 326 20 26
welcome@orfit.com

www.orfit.com



Ref. No. 50257
VERSION 2
LAST UPDATE: 03/08/2022
REVISION DATE: 03/08/2024